BEAUTIFLY

Lipomassage Slim

Body

series

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

 /beautifuly.eu  @beautifuly.eu  beautifuly.eu

 /beautifulyeu  @beautifuly.eu

f1

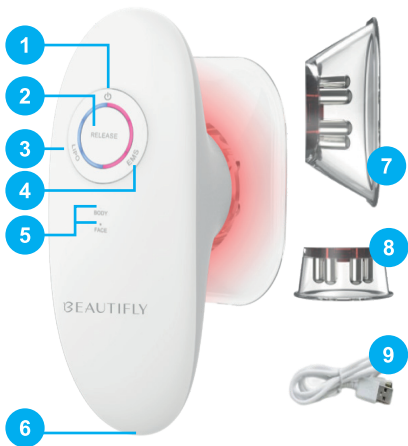f2

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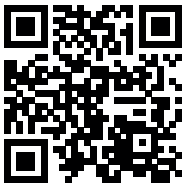- www.beautifly.eu

- www.youtube.com/beautifly.eu

| | |
|-----------------------|-------------|
| Size device: | 145*60*69mm |
| Weight device: | 234g |
| Charging time: | 3-4 hours |
| Working time: | ≥60 minutes |
| Charging method: | Type-c |
| Automatic shutdown: | 10 minutes |
| Red light wavelength: | 630±5nm |

Safety statements

Read the instructions carefully before use. Follow these safety guidelines when using the device.

Instructions for use

Choose the larger or smaller suction cup, depending on the size of the area to be treated. Attach the cup to the device properly locating the lugs and twisting into position.

Device description [page 2 - f1]

1. On/Off button
2. Suction adjustment button
3. LIPO mode button
4. EMS mode button
5. LED indicators
6. Type-C charging input
7. Large suction cup
8. A small suction cup
9. Type-C charging cable

Charging

1. Charge the device by connecting it to a power source using the supplied USB type C cable. While charging, the LEDs will light up sequentially. When the product is fully charged and all six LEDs are lit continuously, disconnect it from the power source.
2. A 2 to 4 hour charge provides up to 1 hour of cordless use.
3. All 6 LEDs will blink 5 times if the battery is low and needs charging.

BODY mode

1. Press the power button to turn on the device, then the device will work in a body mode. LIPO and EMS work on level 1.
2. Press LIPO to increase the suction level, press EMS to increase the EMS level. The suction and EMS operate in 3 levels, indicated by 3 orange LEDs.
Note: In the first use, we suggest starting with level 1 for both the EMS and LIPO modes. Applying a small amount of baby oil or massage oil can help the gliding effect, but avoid using perfumed oils, which can damage the plastic suction cup.
3. Once the settings are set, firmly place the cup on the desired body part and glide across the skin very slowly.
4. For best results, treatment should take 3 to 5 minutes per area every day.
5. When both the LIPO and EMS are on, you should feel both the suction and light electric currents. To stop immediately if either becomes uncomfortable, press the large "Release" button. Adjust the levels to a lower setting before continuing.
6. Press the power button one time to switch to FACE mode or press it twice to turn the device off. The device includes an auto-shutoff after 10 minutes.
7. To remove the cup, twist it counterclockwise.
8. See the 'cleaning' section below for advice on how to change the filter.

FACE mode

1. Press the power button twice and switch to face mode. Only EMS works, so there's no need for a cup.
2. Press EMS to choose the EMS level, 3 levels are indicated by 3 orange LEDs.
3. Once done, press the power button one time to turn off the device.

Tips for Use

- The anti-cellulite massager is especially effective when used regularly and in conjunction with a healthy, balanced diet and regular exercise.
- During usage, keep moving the device on your body. Hold the unit firmly on the body to avoid dropping it. There is no need to push too hard. If the suction power feels uncomfortable or too powerful at any point, push the "RELEASE" button to release the cup from your body.
- Some temporary reddening or bruising of the treated area is not uncommon. This should disappear within a few days.

Cautions

- People with electronic human organs are forbidden from using it.
- Persons with skin surgery within one month or skin sunburn are prohibited from using the device.
- Forbidden to use by a person with damage to the skin, scars, or electromagnetic sensitivity.
- Take off your watch, glasses, and other metal objects before use.
- Avoid shovel head close to the eyeball.
- Don't operate in the same area for a long time.
- Do not use the device during pregnancy.

Cleaning

1. Unplug the device before cleaning.
2. Remove the cup from the main unit.
3. Use tweezers to remove the filter.

Note: We recommend changing the filter every three months or when you start to feel that the suction power is less powerful.

4. After use, the suction cup should be wiped clean with a damp cloth.
5. Never allow water or any other liquids to come into contact with the main unit. Do not immerse in any liquid to clean.
6. Never use abrasive cleaners, brushes, glass or furniture polish, paint thinner etc. to clean the device.

Storage

- Keep it dry and clean.
- Keep it away from high temperatures.
- Keep it out of reach of children.

More information about the company and products [page 2 - f2]

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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