BEAUTIFLY

Derma
Basic
series

Diamond microdermabrasion

User manual

 /beautifly.eu  @beautifly.eu  beautifly.eu

 /beautiflyeu  @beautifly.eu

f1

f2

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- www.beautify.eu

- www.youtube.com/beautifyeu

Power:	0,75W
Working Voltage:	DC 3.7V
Working Amperage:	200mA
Charging voltage:	DC 5V
Dimensions:	4.5x18x5.5cm

Security

Before use, read this manual and follow the recommendations contained therein during use of the device

Device description [page 2 - f1]

1. Head
2. Power button

Landing:

Charge your device for a maximum of 8 hours, making sure it is turned off. Connect via USB cable to a charger or computer and check if the charging indicator (red when charging, turns green when charged) is active. When the battery runs out, the indicator will flash red and green for 10 seconds and then the device will turn off.

How to use:

1. Cleanse and dry your face. This will relax the skin and open the pores. Then place the microcrystalline tip onto the large round cap of the device.
2. Turn on the device and place the tip on your face. Start with your chin and move the device from bottom to top.
3. Place the device on your nose and move from the inside to the outside, thoroughly cleaning each nose lobe.
4. Then place the device on your forehead and move it from the center to the outside in both directions.
5. In the last step, adjust the suction intensity and place the tip on your cheeks. Move the device from the inside to the outside. From the bottom up.
6. After completing the treatment, wash your face with a soothing liquid and use a moisturizer with a filter.

Pore Cleansing:

To open pores and remove blackheads, use peel-off or apply a warm towel to your face for 3-5 minutes. Then select and install the appropriate tip in the device - large or small round - depending on the part and needs of the skin. Aim the tip precisely at the areas to be cleaned by adjusting the suction level. Once done, the toner or mask will close the pores. Move the device slowly, always in the same direction, avoiding oval movements. In case of acne, suck the blackhead and release it immediately.

Wrinkle Removal:

Attach the oval tip to gently massage the area around the eyes, mouth and nose, adjusting the suction level to suit your needs. Tighten the skin and remove wrinkles with low suction intensity. This treatment improves blood circulation and skin elasticity, helping to reduce fine wrinkles. Gently suck and release the skin.

Do not hold in one place for longer than 3 seconds to avoid irritation. In case of irritation, stop using the device for a while.

Precautions and contraindications

Before use, dry your hands and do not rinse or immerse the device under water. Please take care of your facial skin daily when using this device. Do not use other mechanical facial cleansing devices while using this device. Please note that the device has strong suction levels. The heads should be disinfected to avoid the occurrence and spread of viral, bacterial and fungal skin diseases. Contraindications: active purulent acne, rosacea, retinoid therapy, acid treatments, lasers - at least 4 weeks before microdermabrasion, psoriasis, viral, bacterial and fungal skin diseases, vascular skin, tendency to scarring, interrupted skin continuity (fresh wounds), skin cancer or suspicious moles in the area undergoing treatment, fillers in the face - at least 4 weeks of interval before microdermabrasion, otherwise the preparation may move.

Frequency of execution and post-treatment care:

Microdermabrasion is performed every 2 weeks in a series of 6-12 sessions. Depending on the effect we want to achieve, treatment breaks range from 1 to 4 weeks. After the treatment, slight skin peeling may occur for about 3-4 days; the use of additional serum or ampoules is not required. Post-treatment care requires avoiding strong sunlight, chlorinated water, saunas or solariums, and using creams with a high UV filter to protect the skin and increase the effectiveness of the treatment. It is possible to combine it with iontophoresis or needle-free mesotherapy to increase the absorption of active ingredients and achieve better results.

More information about the company and products [page 2 - f2]

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
[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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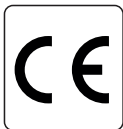[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautifly Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Derma Basic