

BEAUTIFLY

Looky
P R O
series

Fale Radiowe
Rejuvenation Radio
Frequency

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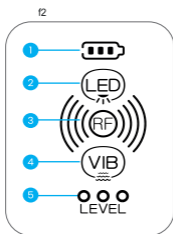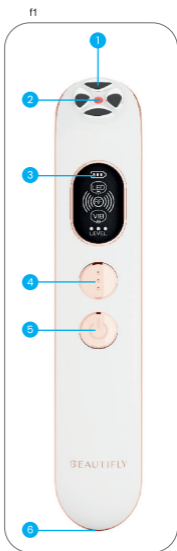- www.beautifly.eu
- www.youtube.com/beautiflyeu

Entry:	DC5V/1A
Battery:	500mAh
Unit weight:	54g
Size:	149.9x28x22mm
Vibration frequency:	12000RPM
RF:	1 MHz
Charging time:	3 hours
Using time:	>30min

Security:

Before using the product, read this user manual and follow the recommendations contained therein when using the device.

Looky PRO is a powerful anti-aging device for the care of the skin around the eyes. Effectively soothes tired and swollen eyes. It stimulates the skin around it, firms it and smoothes it. It has the function of radio waves with a strong lifting effect. High-frequency vibrations of 12,000 RPM combined with red photon light perfectly smooth and care for the thin and sensitive skin around the eyes.

Device description [page 2- f1]:

1. Massage head
2. Red light
3. LCD display
4. Level button
5. Power/Mode button
6. USB charging port

LCD display [page 2- f2]:

1. Battery indicator
2. LED indicator light
3. RF indicator
4. Vibration indicator
5. RF operating levels

Functions:

1. Looky PRO is a device for the care of the eye area. It uses the frequency of radio waves with vibrations, and additionally has a red photon light function. Use it with your eye cream or gel and your skin will become smoother, firmer and brighter.
2. Anti-aging technology using radio waves penetrates the epidermis and heats the dermis. The generated electromagnetic waves and thermal energy exchange each other, and the generated thermal effect promotes the production of collagen and reduces fine wrinkles.
3. Red light treatment increases the skin's capacity for the penetration of active ingredients contained in the care product. Daily use can effectively improve skin elasticity and reduce discoloration around the eyes.
4. High-frequency vibration massage supports blood circulation and better skin absorption of the active ingredients contained in the care product.

Service:

Power button - long press the power button will turn on the device, short press will switch the operating mode.

1. RF + red light + vibration
2. RF + red light
3. RF + vibration
4. RF radio frequency

Intensity level switch button (adjusts RF intensity only)

How to use:

1. Apply an appropriate amount of eye cream to previously cleansed skin. Spread it evenly.
2. Press the power button to turn on the device, short press the power button to switch the working mode, and short press the power button to adjust the RF intensity.
3. Use the massage head as shown for massage in the drawing.
4. After use, wipe the massage head with a dry or paper towel.

Storage:

1. Eye care products should be distributed evenly before each treatment to avoid excessive accumulation of eye care products in the massage head.
2. After each treatment, wipe the product with a dry towel.
3. Do not wash the device under water or submerge it. The device is not waterproof.
4. Do not use corrosive agents to clean the device.
5. Store it in a cool, dry place.
6. Do not use the device if it is damaged.
7. Do not press the device against the skin with great force.

Charging method:

1. At low power level, turn off the device and connect to charge.
2. During the charging process, the battery symbol on the display flashes.
3. Charge your device regularly to avoid battery aging.
4. Once charging is complete, disconnect the charger from the device.

Precautions:

1. Follow the recommended operating instructions. If your eyes have inflammation such as conjunctivitis, you should not use this device until it is completely healed. Do not use this product while wearing contact lenses.
2. Remove all metal elements from the facial area before using the massager, such as earrings and necklaces.
3. Do not use the Looky PRO device while charging.
4. Do not repair, disassemble or modify the device yourself. Any irregularities should be assessed by the service.
5. Please do not drop or hit the device to avoid it malfunctioning.

6. When the power cord is damaged, do not use it to avoid short circuit or fire.

7. The device should not be used on people:

- with infectious diseases
- with heart disease
- with high fever
- with malignant tumor
- with acute inflammation
- with a pacemaker
- with skin diseases
- pregnant/breastfeeding
- with excessive sunburn

More information about the company and products [page 2 - f3]:

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[EN] Rules for handling electrical and electronic equipment: Electrical and electronic equipment as well as batteries and accumulators must not be combined with other household waste. This is indicated by the crossed-out wheeled bin symbol on the packaging. Dispose of waste electrical and electronic equipment as well as batteries and accumulators at your local separate collection point or to your dealer. Electrical and electronic equipment and batteries and accumulators may contain hazardous substances, mixtures and components that are harmful to the environment and human health. Incorrect handling or damage to them may lead to damage to health or environmental pollution when the equipment is further disposed of or recycled. The device was placed on the market after August 13, 2005

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[EN] Summary declaration of conformity: The manufacturer hereby declares under its sole responsibility that the above product complies with Directive 2014/30/EU of the European Parliament and of the Council of 26 February 2014 on the harmonisation of the laws of the Member States relating to electromagnetic compatibility and Directive 2011/65/EU of the European Parliament and of the Council of 8 June 2011 on the restriction of the use of certain hazardous substances in electrical and electronic equipment, together with Commission Delegated Directive (EU) 2015/863 of 31 March 2015 amending Annex II to Directive 2011/65/EU of the European Parliament and of the Council as regards the list of restricted substances.

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Importer: Beautifly Sp. z o.o.

Plac Bankowy 2, 00-095 Warszawa, Polska

Model: Looky PRO